



REHABILITATION MEDICAL COMPANY



Experience The Technology of **No.1** Electroceutical **NMS** Therapy Provided by **REMED**

Mission


Improve the field of electroceutical worldwide as a No.1 leading company with outstanding technology and solutions.

Vision

Offer Brain to Body care for people everywhere.

Values

Passion, Innovation, Technology, Trust

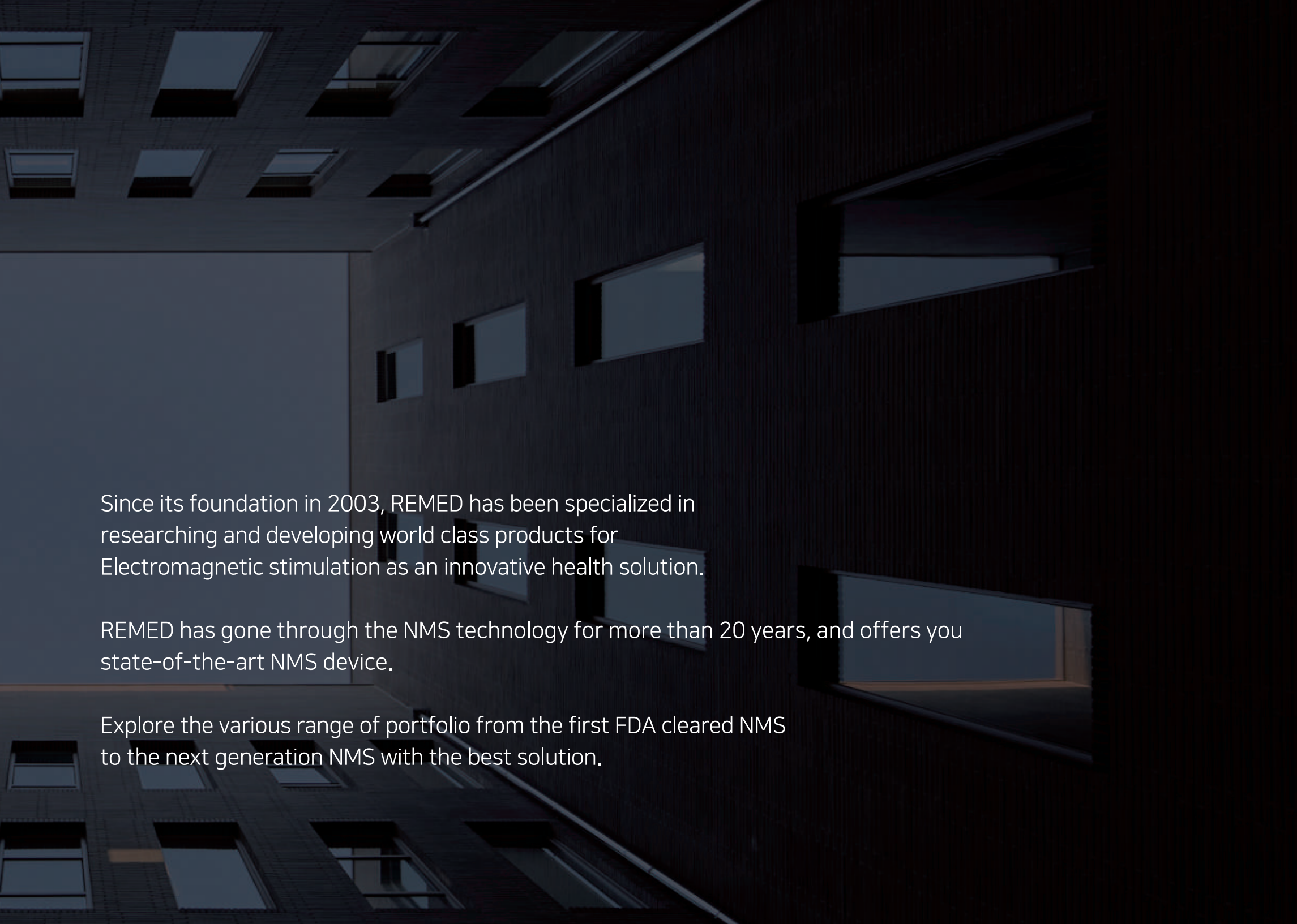

Since its foundation in 2003, REMED has been specialized in researching and developing world class products for Electromagnetic stimulation as an innovative health solution.

REMED has gone through the NMS technology for more than 20 years, and offers you state-of-the-art NMS device.

Explore the various range of portfolio from the first FDA cleared NMS to the next generation NMS with the best solution.

Exporting to over **50** countries

- 📍 REMED
- 📍 REMED Subsidiary
- Distributors


REMED NMS

system has been installed more than 90% in domestic market since 2007. It has enlarged its market worldwide and has successfully been introduced to the world.







WHAT IS PAIN MANAGEMENT SOLUTION?

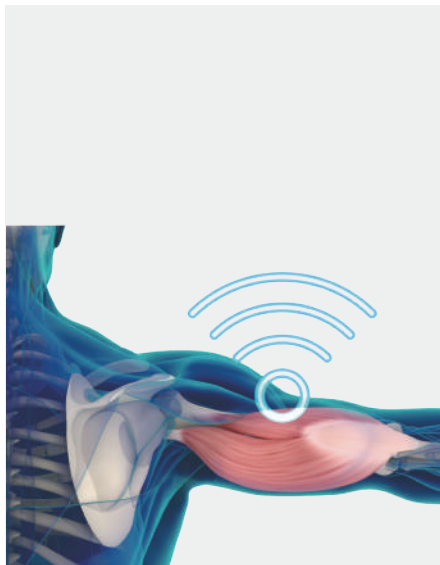

FSWT

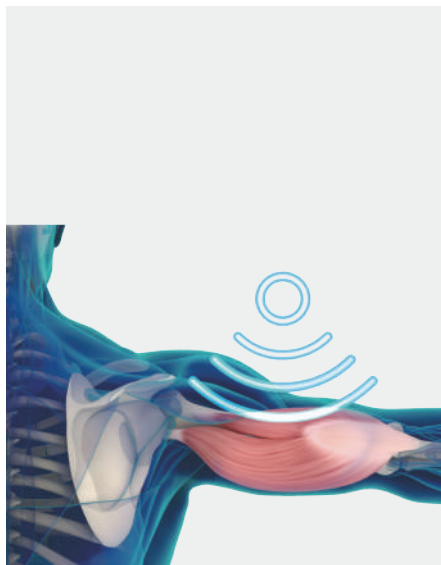

RSWT

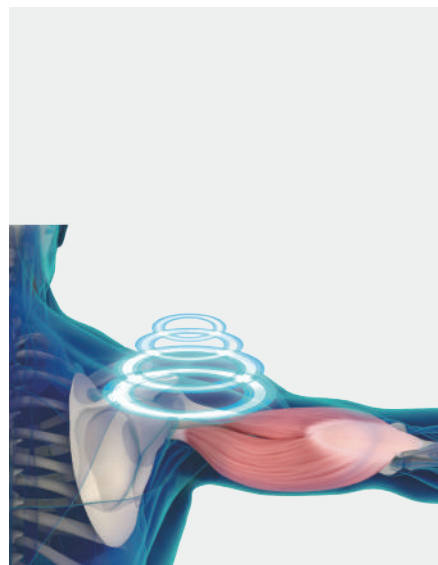

NMS

PMS stands for Pain Management Solution provided by REMED, and it is a way of relieving pain with three different types of energy sources. These three energy sources are FSWT (Focused Shock Wave Therapy), RSWT (Radial Shock Wave Therapy), and NMS (Neuro magnetic stimulation), which can be applied to patients suffering from acute and chronic pains in the musculoskeletal system. These three different energy sources generate various clinical effects on the human body, and they eventually lead to tissue regeneration, which means these treatments restore cells, blood vessels, tendons, and ligaments.



FSWT treatment
less than 5 mins



RSWT treatment
less than 5 mins



NMS treatment
less than 15 mins



**Regenerative Pain
Treatment**



EXTRACORPOREAL SHOCK WAVE THERAPY SOLUTION



SALUS-RSWT

1CH Radial

SALUS-ESWT-EXPERT

1CH Radial &
1CH Focused

SALUS-FSWT-DUET

2CH Radial &
1CH Focused

OPTIMUS-PRO

1CH Radial &
1CH Electro Magnetic

WHAT IS EXTRACORPOREAL SHOCK WAVE THERAPY?

Extracorporeal Shock Wave Therapy (ESWT) is a non-invasive treatment using acoustic shock waves to enhance healing and reduce pain in musculoskeletal conditions. It increases blood flow, breaks down calcifications, and promotes new blood vessel formation. ESWT is used for conditions like plantar fasciitis, tennis elbow, shoulder tendinopathy, and other chronic soft tissue injuries.



Energy Level	1~6 Bar, 0.1 Bar Step	0.5~6 Bar, 0.1 Bar Step Max 1.0mJ/mm ²	1~6 Bar, 0.1 Bar Step Max 1.0mJ/mm ²	3 Tesla / 1~6 Bar, 0.1 Bar Step
Frequency	1~21Hz, 1 Hz Step	1~21Hz / 1~12 Hz, 1 Hz Step	1~21Hz / 1~12 Hz, 1 Hz Step	1~100Hz / 1~21 Hz, 1 Hz Step
Channel	1 CH / 2 CH (Option)	2 CH (RSWT + FSWT)	2 CH (RSWT + FSWT)	2 CH (Electromagnetic + RSWT)
Treatment mode	M1~M20	M1~M20 / 1~20 Protocols	M1~M20 / 1~20 Protocols	A1~A20, M1~M20, U1~U20 / M1~M20
Applicator	6, 10, 15, (20 & 35mm, Option)	6, 10, 15, (20 & 35mm, Option) 5, 10, 15, 20, 25, 30, 35, 40 mm 0, 5, 10, 15, 20 (Linear, Option)	6, 10, 15, 20 & 35mm (Option) 5, 10, 15, 20, 25, 30, 35, 40 mm 0, 5, 10, 15, 20 (Linear, Option)	1 Transducer 6, 10, 15, 20 & 35mm (Option)
Interface	8" LCD Touch Screen, Jog Shuttle	10.2" LCD Touch Screen, Jog Shuttle	10.2" LCD Touch Screen, Jog Shuttle	8" LCD Touch Screen, Jog Shuttle
Dimension	47 X 50 X 99 cm	51.6 X 48.1 X 114.5 cm	66 X 64 X 104 cm	47 X 60 X 110 cm

ELECTROMAGNETIC STIMULATION SOLUTION



SALUS-TALENT
Basic Line



SALUS-TALENT-A
Standard Line



SALUS-TALENT-Pro
Professional Line



PERFECT
Premium Line

WHAT IS ELECTROMAGNETIC STIMULATION?

NMS Line generate magnetic field of 3 Tesla which is about 600 times stronger than the normal bar magnet and penetrates deep inside of the body. This strong electromagnetic power has made physiotherapists provide highly effective and fundamental treatment for acute and chronic pains.



Model
Type
Channel
Magnetic Strength
Magnetic Frequency
Treatment Mode
Interface
Dimension

SALUS-TALENT
 ELECTROMAGNETIC
 1CH
 2.0 Tesla
 1~50Hz
 Auto : 4, Manual : 4
 6 Buttons, 1 Jog dial
 40 X 54 X 119 cm


SALUS-TALENT-A
 ELECTROMAGNETIC
 1CH
 3 Tesla
 1~100Hz
 Auto : 9, Manual : 9
 4 Buttons, 1 Jog dial
 46 X 53 X 98 cm

SALUS-TALENT-Pro
 ELECTROMAGNETIC
 2CH
 3 Tesla
 1~100Hz
 Auto : 20, Manual : 20, User : 20
 8" LCD Touch Screen, 1 Jog dial
 55 X 53 X 100 cm

PERFECT
 ELECTROMAGNETIC
 2CH (SIMULTANEOUS)
 3 Tesla
 1~200Hz
 Specialized area mode
 12.1" LCD Touch Screen
 54 x 66 x 112 cm

REMED APPLICATORS

NMS

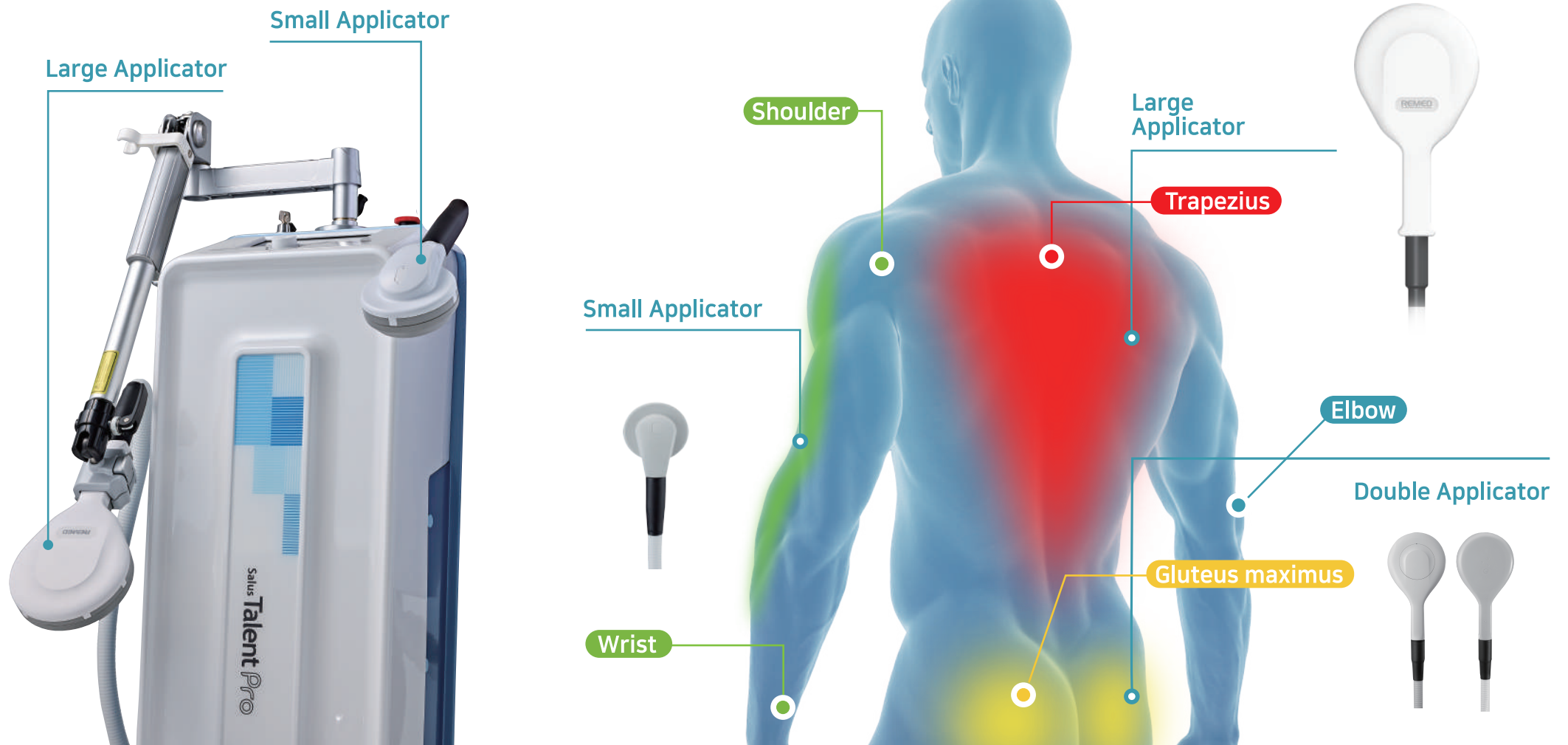

ESWT RSWT



S6.

Acupuncture points



S10.

Trigger points, deep-lying enthesopathy



S15.

Any type of enthesopathy



S20.

Muscle tissue and connective tissue



S35.

Muscle tissue and connective tissue, large muscles



T15.

Deep-lying pathologies, chronic conditions, high-energy application



T20.

Muscle tissue and connective tissue, high-energy application



REMED APPLICATORS

ESWT FSWT



R-PAD-5
Depth 0.5cm



R-PAD-10
Depth 1cm



R-PAD-15
Depth 1.5cm



R-PAD-20
Depth 2cm



R-PAD-25
Depth 2.5cm



R-PAD-30
Depth 3cm



R-PAD-35
Depth 3.5cm



R-PAD-40
Depth 4cm



L-PAD-0
Depth 0cm



L-PAD-5
Depth 0.5cm




L-PAD-15
Depth 1.5cm



L-PAD-10
Depth 1cm




L-PAD-20
Depth 2cm

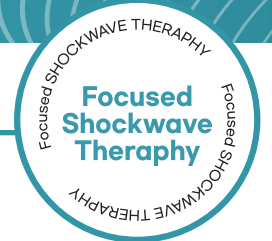

Joint | Tendinitis | Ligament

Superficial | Tendinitis

INDICATION

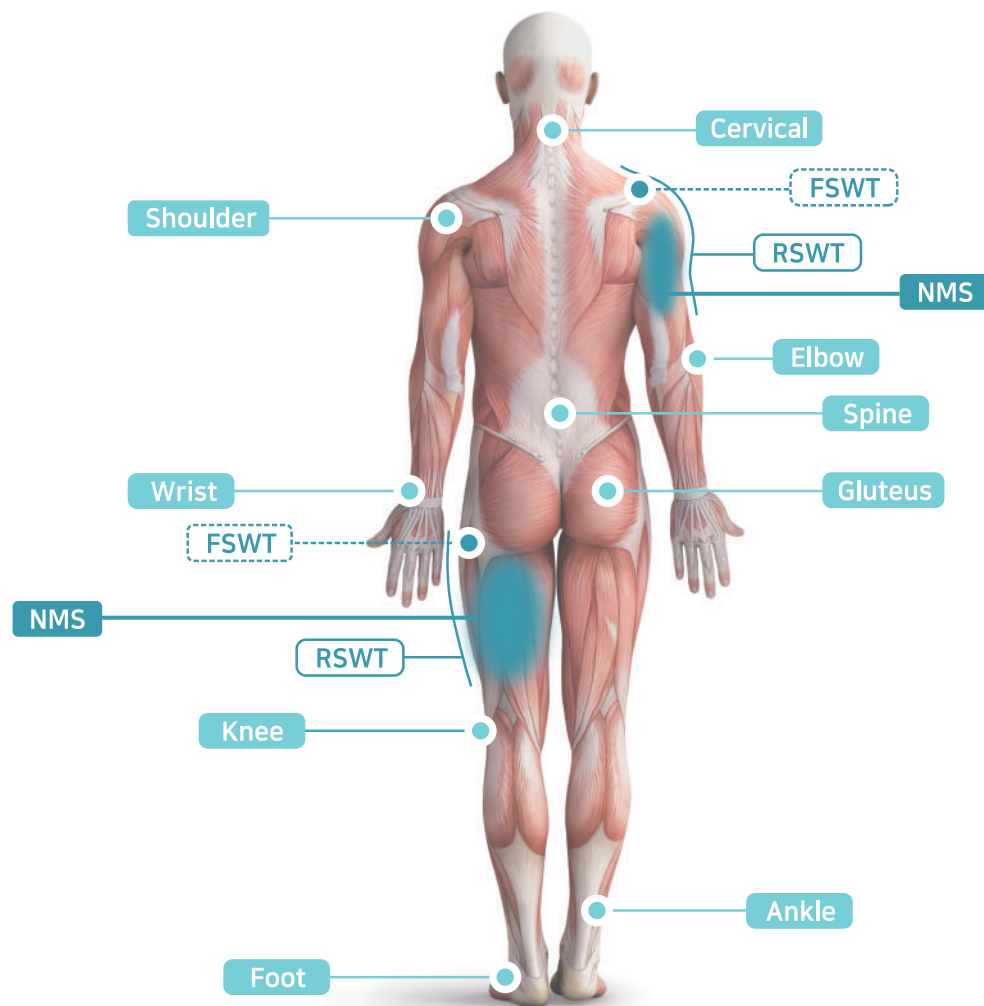

SALUS-TALENT SERIES

Musculoskeletal Disorders

Degenerative Arthritis, Rheumatoid Arthritis, Cervical Pain, Sport Injuries, Frozen Shoulder, Muscle Relaxation, Pain Control After Car Accident

Spinal Disorders

Acute/Chronic Lumbago, Hip Gout, Sciatica, Spinal Bifida, Spondylitis, Chronic Back Pain, Spinal Disc Disease, IMS Treatment, TPI Treatment

SALUS-RSWT, OPTIMUS-PRO

Painful Shoulder, Calcified Shoulder, Frozen Shoulder, Plantar Fasciitis, Tennis & Golfer's Elbow, Achilles Tendon, Delayed Bone Fracture Recovery, Pain for knee, Back & Neck Degenerative Arthritis, Muscle Relaxation, Spinal Stenosis, Calcific Tendinitis, Pelvic Pain, Joint Pain, Acute or Chronic Pain

SALUS-FSWT, SALUS-ESWT-EXPERT

Thumb Pain, Wrist, Epicondylitis(Ulnar), Epicondylitis(Radial), Patellar Tendinitis, Ankle, Achilles Tendinitis, Plantar Fasciitis, Cervical & Shoulder Pain, Calcific Tendinitis, Supraspinatus Tendon, Myofascial Pain, HIP & Lumbar & Trochanter Pain, Thigh Treatment

NMS

: Neuro Magnetic Stimulation

NMS



REMED Pain Treatment In 4 steps



Pectoralis minor muscle with NMS



Gluteus maximus muscle with NMS



Low Back Pain with NMS



Upper trapezius muscle with NMS

ESWT

: Extracorporeal Shock Wave Therapy

ESWT



REMED Pain Treatment In 4 steps

01



Treatment of Rotator cuff with FSWT(Focused)

02



Treatment of Tennis elbow with FSWT(Linear)



Low-Back Pain with RSWT



Soleus muscle with RSWT

WHERE THE TECHNOLOGY & INNOVATION MEET

Excellence in NMS technology over 20 years, REMED is developing together with customers based on the deep trust as we pursue the shared value and virtuous cycle of the world.



2003 COMPANY FOUND
REMED was established in the year of 2003.

2008 THE 1st NMS SYSTEM
The first NMS system launched and developed by REMED

2014 CE MARK CLEARED
The first Asian CE mark cleared TMS & NMS system.



2006 OIL COOLING SYSTEM PATENT
Patent for the first Oil Cooling System for effective coil cooling.



2013 WORLD FDA
Cleared FDA from Korea, China, Russia, Thailand, etc.



2015 WORLD CLASS PRODUCT
Awarded as the next promising World Class Product in Korea.

2016 JFDA CLEARED
Japanese FDA cleared as
the First Asian TMS & NMS system.

KOSDAQ

2019 KOSDAQ LISTED
Listed on KOSDAQ
stock market.

2021 US FDA CLEARED
Salus-Talent-Pro
FDA approved



20
ANNIVERSARY

2023 20 YEARS ANNIVERSARY
REMEDI has celebrated 20 years of foundation as leading TMS & NMS company.



2016 ISO13485 CERTIFIED



2020 EXPORT of
10 MILLION DOLLARS
Awarded by Ministry of Industry
as the next leading company.

REMEDI
JAPAN

2022 REMEDI JAPAN
The first Subsidiary, foreign branch
REMEDI Japan established.



2024 DEVELOPED THE NEXT
GENERATION TMS
Finally the newest next TMS system
developed and launched.



Tel +82 31 606 5761 | Fax +82 31 606 5757 | E mail kdy1741@remed.kr | Website www.remed.kr

Sales Office / R&D Center 21-7, Wiryeseoil-ro 1-gil, Sujeong-gu, Seongnam-si, Gyeonggi-do, Republic of Korea 13647

Ver.03

